

# Youmount®

CATHETER MOUNT

**CORRUGATED  
AND EXTENSIBLE  
CATHETER  
MOUNT**



The **CATHETER MOUNT YOUMOUNT®** is a corrugated tube, indicated for the transport of anesthetic gases or oxygen to patients undergoing mechanical ventilation.

The **CATHETER MOUNT YOUMOUNT®** is connected between the patient (endotracheal tube or tracheostomy) and the breathing circuit.




### **Biconnect Cap (dual function):**

It can be used at the time of intubation and at the time of filter disposal.

## Cateter Mount

- Corrugated and extensible tube
- ISO standard conical connections, 22/15 mm
- Dual-purpose tube cap

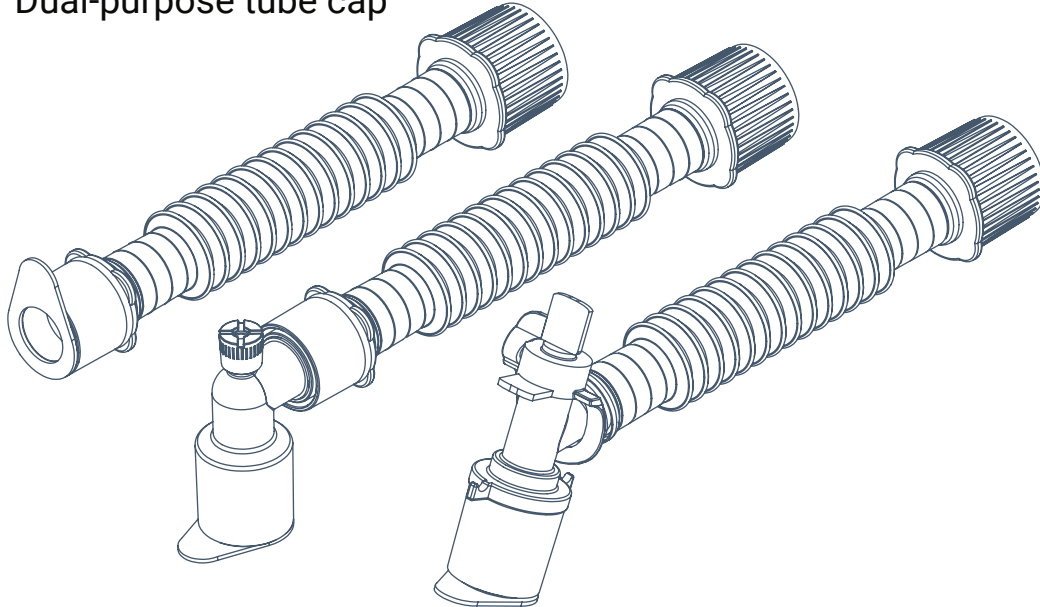

The **CATHETER MOUNT YOU MOUNT®** should be used according to its presentation:

- **YMO-R (STRAIGHT)**: In connection with the breathing filter in tracheostomized patients or with endotracheal intubation;
- **YMO-A (ANGLED)**: In connection with the breathing filter in tracheostomized patients or with endotracheal intubation. This model has connection with Luer-Lock® cap for capnography.
- **YMO-AG (ANGLED ROTATING)**: In connection with the breathing filter in tracheostomized patients or with endotracheal intubation. This model has a angled rotating connector that moves with the circuit, and a silicone top cap that allows the insertion of aspiration catheters or bronchoscope.

### RETO

- Tamanho: 17 cm
- Conexões: 15F/22M \* 22F

### ANGULADO

- Tamanho: 23 cm
- Conexões: 15M 90° \* 22F \* capnografia

### ANGLED ROTATING

- Size 19 cm
- Connections: 15F/22M \* 22F \* Silicon Cap for Aspiration Catheter or Bronchoscope

| Product/Description     | Code   | Anvisa      |
|-------------------------|--------|-------------|
| Straight Catheter Mount | YMO-R  | 81484760013 |
| Angled Catheter Mount   | YMO-A  | 81484760013 |
| Rotating Catheter Mount | YMO-AG | 81484760013 |